


YSI 5000 Series Products

Advanced Dissolved Oxygen Instruments for the Lab

Ease of use, data integrity, expanded range, enhanced digital capabilities – the YSI 5000 Series is engineered to meet the need of highly accurate, reliable dissolved oxygen measurements, primarily to determine BOD, Oxygen Uptake Rate (OUR) and Specific Oxygen Uptake Rate (SOUR).

YSI 5000 and 5100 Instruments

The YSI 5000 features menu-driven operation and a large liquid crystal display. Internal memory for 100 data sets with date/time stamp eliminates the need for manual entry of readings. An RS-232 interface allows upload to spreadsheets. Multitasking is enhanced with the auto-stabilization capability, which enables you to predetermine a stable reading and be alerted to that condition. The full-featured YSI 5100 instrument offers even greater performance and flexibility including built-in software for OUR/SOUR measurement.

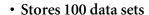

- RS-232 interface for direct link to your computer
 - Large graphic display and keyboard profile
 - Bar code capable*
 - Built-in barometer for autocalibration*
 - On-board software for OUR/SOUR measurement*

Pure Data for a Healthy Planet."

5100 pictured

Saves time

Automates BOD calculations

Reduces errors

Insures data integrity

Conforms to Standard Methods

YSI 5120 BOD Analyst™ software for Windows ®

- Eliminates need for initial measurements of each sample dilution.
- Conforms to Standard Methods for seed depletion, minimum final DO, and minimum depletion levels.
- Built-in safeguards include notification if seeded samples are present without the corresponding seed control samples.
- Sample Groups allow assignment of a single name to a number of samples, reducing setup time.
- * Only available on the YSI 5100

Y S I Environmental



Pure Data for a **Healthy**

To order or for more information, contact **YSI Environmental.**

800 897-4151

www.YSI.com

YSI Environmental 937 767 7241 Fax 937 767 9353 support@YSI.com

Endeco/YSI 508 748 0366 Fax 508 748 2543 support@YSI.com

YSI Environmental European Support Centre 44 1730 710 615 Fax 44 1730 710 614 europe@YSI.com

YSI (Hong Kong) Limited 852 2891 8154 Fax 852 2834 0034 hongkong@YSI.com

YSI/Nanotech (Japan) 81 44 222 0009 Fax 81 44 222 1102 nanotech@YSI.com

YSI (Qingdao) Limited 86 532 389 6648 Fax 86 532 389 6647 china@YSI.com



BOD Analyst, Who's Minding the Planet? and Pure Data for a Healthy Planet are trademarks of YSI Incorporated. Windows is a registered trademark of Microsoft Corpora

Printed in USA 0303 W21-02 CE

> YSI incorporated Who's Minding the Planet?™

For more information on YSI's full line of products and services, as well as application notes, FAQs, troubleshooting, and more, check our information-packed website: www.YSl.com

YSI 5000 Features Benefits • Ease of use • Menu-driven operation · Compatible with all YSI probes · Less investment needed · Easy-to-read from 10 feet · Large liquid crystal display · Internal memory for 100 data sets · No manual transcribing • Compatible with YSI BOD Analyst™ software · Auto BOD calculations; no errors · RS-232 interface · Compatible with PC lab information systems • Probe power through instrument · Fewer cables in workspace · Control functions from keyboard Computer interface control · Audio & tactile keyboard response · Keystroke confirmation reduces errors · Autostabilization • Alerts to predetermined stable condition • Ensures upgrade without downtime · User-upgradable internal software Date-stamped stored data · Real-time clock

TSI STUU (provides these additional features & benefits)				
Features	Benefits			
Built-in OUR/SOUR software	 Eases use with USEPA 503 regulations & for operational control 			
Internal barometer	 No external barometer or look-up charts required 			
Automated calibration	Reduces manual steps; saves time			
Barcode compatibility	Automatic bottle numbering			
Port for computer keyboard	Quick data entry from keyboard			

YSI 5000 & 5100 Specifications

Accuracy	mg/L	±0.1% + 1 lsd	Resolution	mg/L	0.01
	% air	$\pm 0.1\% + 1 \text{ lsd}$		% air	0.1%
	temperature	±0.1° C		temperature	0.1° C
Range	mg/L	0-60	Power	battery & AC	
	% air	0-600%			
	temperature	-5 to +50° C			

YSI 5120 BOD Analyst™ for Windows®

Features

- Normalize to a 5, 7, 20, or 28 day BOD, when final readings are early or late.
- Read data directly from the YSI 5000 or 5100 Dissolved Oxygen Instrument or enter it from the keyboard in spreadsheet format.
- · Make graphs more meaningful by plotting statistics that are of interest.
- Eliminate errors by requiring bottle number entry via barcode scanner only (5100 only). Bottle numbers can also be entered via the YSI 5000 or 5100 keypad or a PC keyboard.

How to order

- YSI 5000 Laboratory Dissolved Oxygen (115V)
- YSI 5100 DO Lab Meter with Built-in OUR/SOUR Software (115V)
- · YSI 5015 Barcode Scanner
- YSI 059160 Barcode Labels
- YSI 5120 BOD Analyst[™] for Windows[®]
- To order probes that meet the needs of your application, see your YSI representative.